

breeze™
Easy Digital



unitron™
hearing 

DIGITAL ADVANTAGES MADE EASY

The benefits of a digital instrument, the simplicity of a conventional fitting

Offering the advantages of digital hearing instruments to everyone is now easier and more affordable than ever. Breeze is the easy choice for those seeking the benefits of digital hearing instruments, who have previously been constrained by lack of access to computer technology for fitting. Breeze combines the benefits of digital sound quality with the ease of conventional, trimmer-controlled fittings.

The Benefits of Digital Sound Processing

- ▶ Clear, comfortable sound
- ▶ Reduced background noise
- ▶ A more precise fitting

The Simplicity of Conventional Fittings

- ▶ No computer, cables or software required
- ▶ Simple, effective trimmer-controlled fitting
- ▶ A personalized fit anywhere

Breeze makes digital easy: easy to fit, easy to afford, easy to use.






THERE'S A BREEZE FOR EVERYONE

Choice of processing strategies provides the right fit for everyone

Breeze addresses a wide range of hearing losses, from mild to severe, allowing you to satisfy a wider range of clients simply and quickly. When ordering, select one of three pre-set processing strategies – WDRC, AGCo or Linear – depending on your clients' individual hearing losses and preferences.


WDRC

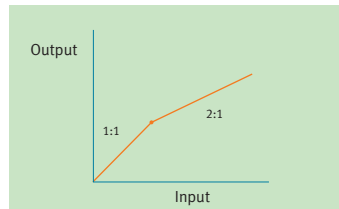

Wide Dynamic Range Compression
For mild to moderately severe hearing losses.

Breeze WDRC allows for restoration of normal loudness growth, permitting wearers to hear soft sounds they may have previously missed.


AGCo

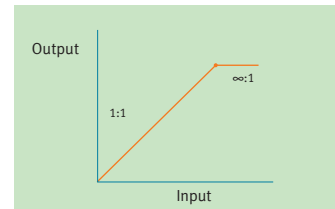

Output Compression Limiting
For mild to severe hearing losses.

Breeze AGCo offers the advantages of compression for louder inputs, and is appropriate for previous linear wearers or more severe hearing losses.

Linear

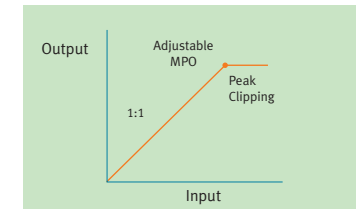

Linear Limiting
For mild to severe hearing losses.

Breeze Linear is an excellent way to introduce analog linear wearers to digital sound.

BREEZE FITS ANYWHERE

Simple trimmer controls yield an individualized fit on the spot


Lack of computer access for fitting is no longer a barrier to offering personalized digital hearing instruments to your clients. Achieve a quick, accurate fit with Breeze's simple, trimmer-based controls. As you rotate the trimmer controls, your clients receive an audible signal of the adjustments. This allows them to provide immediate input regarding their preferences, yielding an individualized fit.

Longer battery life with power management system

Breeze helps ensure a longer battery life thanks to a special power management system that uses power efficiently. In addition, when the battery becomes low, Breeze conveniently emits a series of discreet beeps, giving clients advance notification.

Clean sound with quiet mode expansion

Hearing aid wearers often complain of perceived noise in the instrument, whether from microphone noise or the amplification of low-level sounds in the environment, such as the hum of a fan. Breeze addresses this concern with quiet mode expansion, which reduces gain for very soft inputs, yet preserves moderately soft inputs such as speech, allowing your clients to enjoy more natural, clean sound.

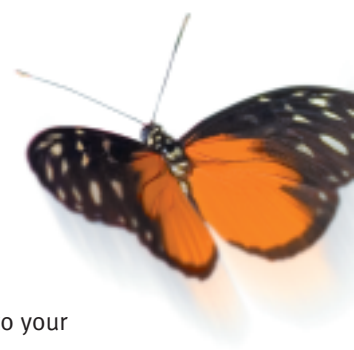



THE EASY CHOICE

Making the move to digital has never been easier.

Breeze delivers the advantages of digital sound processing to everyone on your client list.

- ▶ Digital sound quality
- ▶ Simple, trimmer-controlled fittings
- ▶ Individualized fit
- ▶ Choice of processing schemes address wide range of losses
- ▶ Convenient features

About Unitron Hearing

Unitron Hearing is committed to bettering the lives of people with hearing loss by developing high-quality hearing solutions that incorporate practical features to solve the everyday problems and concerns they have with hearing loss and hearing aids. This commitment also benefits you, the hearing healthcare professional, with practical fitting options, support when you need it and the ability to offer solutions that address what matters most to your clients.



www.unitronhearing.com